




# Surgical Mask Type II

EN14683

TYPE II

0%  
latex

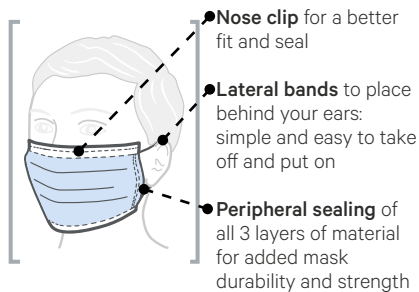

## RESPIRATORY PROTECTION DESCRIPTION AND COMPOSITION:

**Type II surgical mask** (bacterial protection) with lateral bands for easy positioning behind your ears and for a secure fit, making it easy to take off and put on. It comes with a **nose adapter** to increase the mask's fit.

**It consists of 3 layers of material: the outer layers are made of non-woven polypropylene and an intermediate filtering layer.**

Its **≥99% filtering efficiency is based on the BFE test** (Bacterial Filtration Efficiency), exceeding 98% which is the minimum required by the EN14683 standard for Type II masks.

These are **single-use surgical masks**, so they must be disposed of once used.

These masks **can be used in combination with other protective equipment** (glasses, visors, caps, etc.) for overall user protection.

The mask measures 17,5cm x 9,5cm (±0,5cm).

**Latex, Metal and Silicone-free**



**Non-sterile**



**Do not reuse:**  
Single-use mask



**TYPE II** surgical masks  
that offer a bacterial  
filtration of: ≥99% (BFE)


Ref.	Product	BFE
914220	Box of 50 surgical masks	≥99%

# RESPIRATORY PROTECTION

MARKING INFORMATION				
Standard and Certification	EN14683 Surgical masks			
Requirements Standard Minimums	Bacterial Filtering Efficiency (%)	Breathability or Differential Pressure (Pa/cm²)	Fluid Resistance (kPa)	Microbial Cleansing (cfu/g)
	≥ 98	< 29.4	Not required	≤ 30

COMPONENTS			
Latex	Not present	Materials:	Polypropylene 23 gr. + Meltblown 25 gr. + Polypropylene 25 gr  Nasal bridge: Plastic+Iron core
Metal	Not present		
Silicone	Not present*		

\* Although the mask does not contain silicone, we recommend the product be tested if it is used in areas where components may react to this material's presence.

OTHER FEATURES					
Uses	Surgical masks designed to limit the spread of infectious agents from staff to patients during surgical procedures in medical services. It is not a PPE, nor does it provide protection against particulates. Non-sterile. Do not reuse.				
Conservation Storage Expiry	Store in a storage case in a cool, well ventilated place, avoiding dirt and dust. ☀️ Keep away from heat sources. ☔ Protect from rain and humidity.				
Instructions Use How to use	This equipment offers very high bacterial protection, whenever used correctly. To do so, follow these instructions: <ul style="list-style-type: none"> <li>- Wash your hands with alcohol-based disinfectants or soap and water for at least 40-60 seconds before putting on the mask.</li> <li>- Cover your nose and mouth with the mask and make sure you place the mask as close to your face as possible.</li> <li>- Avoid touching the mask during use.</li> <li>- Replace the mask with a new one when it becomes wet.</li> <li>- Discard the mask after usage.</li> <li>- Wash hands again with alcohol-based disinfectants or soap and water for at least 40-60 seconds after use and after disposal of the mask.</li> </ul>				
Presentation	50 units per box. 40 box carton. Carton size: L-52cm x W-38cm x H-30cm 				
Bar code	<table> <tr> <th>BOX BAR CODE</th><th>CARTON BAR CODE</th></tr> <tr> <td>8423173893378</td><td>18423173893375</td></tr> </table>	BOX BAR CODE	CARTON BAR CODE	8423173893378	18423173893375
BOX BAR CODE	CARTON BAR CODE				
8423173893378	18423173893375				